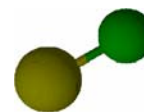


S2355NN **Gold Line****Typical analysis**

Assay (Argentimetric)	98%
Water solubility	0.5%

Packaging

10g black plastic jar

**Physical properties, composition and data**

Chemical formula	AgBr
Chemical group	Bromides
Atomic weight	187.77
Appearance	Pale yellow crystals or yellowish powder
Solubility	Slightly soluble in water, soluble in dilute ammonia
Other information	Darkens on exposure to light

Laboratory preparation, applications and practices

Laboratory preparation Silver nitrate is dissolved in water and a solution of alkali bromide added slowly. The precipitated silver bromide is washed repeatedly with hot water; the operation must be carried out in a darkroom under a ruby red light

Usage Analytical reagent

Storage and handling information

Storage	Protect from light
Disposal	13

Transport regulations

Tariff code	2843.29.00
-------------	------------

